

Neonatal Suction Regulators

In Suction Management, We Are
The Calm Inside the Storm.

Features Specifically Designed for the NICU
Nurse and Patient

- Upper suction limit meets current clinical practice guidelines for infants
- High flow at low pressure for safe but efficient secretion removal
- Full 12-year factory warranty
- Self-Clearing technology™
- Accurate at low suction levels
- Sterilizable - VSS270



Linear Gauge


- As suction increases, the piston goes down indicating the level of suction with a 180 field of view
- Color coding readily shows common levels of suction

Mode Selector Valve

- Easily switch between operating modes OFF-INTERMIT-CONT

Adjusting Knob

- Turn COUNTER-CLOCKWISE to increase the suction setting
- Turn CLOCKWISE to decrease suction setting
- Similar function to flowmeters

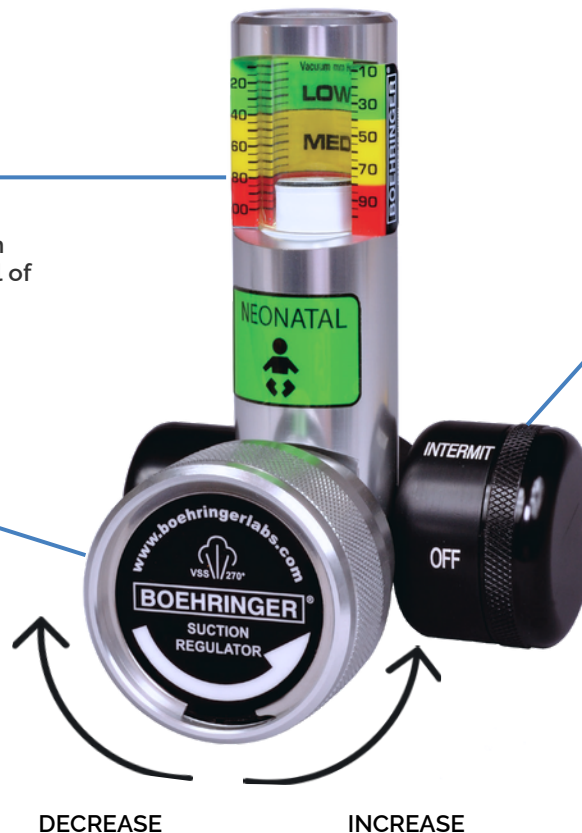

Description	Part Number	Ordering
Neonatal Suction Regulator	3814	1 unit